


Quality-Water

SPECIALISTS

Sampler Audit Booklet

Sampler Audit Booklet

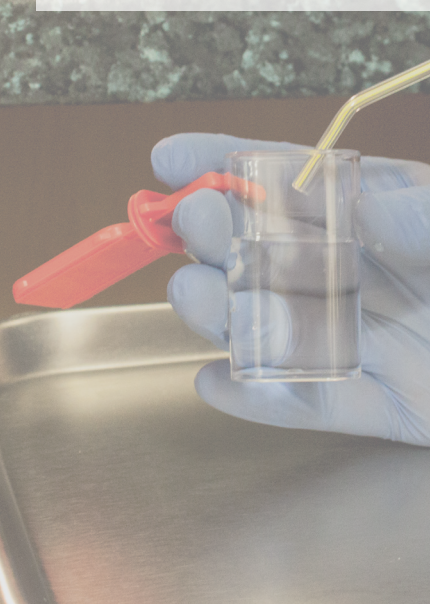

Aura Infection Control Ltd
14b Redwell Court
Harmire Enterprise Park
Barnard Castle
County Durham
DL12 8BN

VAT Reg: GB 181 4500 29
Company Number: 08814266

01833 630393
www.aiconline.co.uk

Auditing > Reading sampler results

To understand how to read the results of the Samplers follow the flow diagram shown here. If you are not sure about your result please don't hesitate to call our clinical support team

who will be happy to cast an eye over your result and give you some advice.

Step 1 Check Sampler after 7 days

It is important to check after exactly 7 days for an accurate result. Use a mag lamp for accurate inspection. Remember you're looking for distinct colonies or a complete overgrowth where all colonies have formed together.
Is the Sampler 100% Clear, with a matt black Sampler surface?

No

The Sampler is not clear or has a shiny surface

Proceed to step 2.

Yes

The Sampler is 100% clear

Record 0 C.F.U. on your audit sheet. Audit complete.

Step 2 Count and record the colonies

Now you need to count the colonies (dots). Use a mag lamp. Record the number of colonies as a C.F.U. on your audit sheet. If there are 'too many to count' or the sampler is slimy record, TMTC. Proceed to step 3.

Step 3 Mark as a pass or fail

Now determine whether the audit has passed or failed using the C.F.U.

Pass

C.F.U. is less than 99

Mark audit as a pass. Audit complete.

Fail

C.F.U. is more than 99

Mark audit as a fail. Audit complete.

Fail

Total overgrowth or slimy


Mark audit as a fail. Audit complete.

Auditing > Sampler examples

E.g. 1 C.F.U.

A red plastic sampler with a black grid. A single white dot is visible on the grid.

E.g. 41 C.F.U.

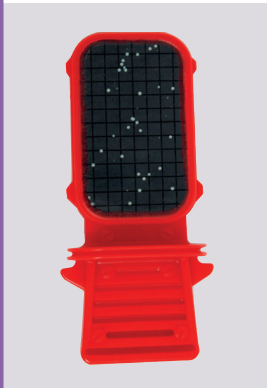A red plastic sampler with a black grid. Approximately 41 small white dots are visible on the grid.

E.g. 100+ countable C.F.U.

A red plastic sampler with a black grid. The grid is covered with many small white dots, representing more than 100 colonies.

E.g. TMTC Total overgrowth

A red plastic sampler with a black grid. The grid is completely covered with a thick, white, fuzzy growth, representing total overgrowth.

Audit conditions > Caution Samplers

A caution refers to a Sampler that has between 80-99 colonies. Although this is a pass, we recommend you look into why the count is high, especially if the count was previously lower.

There can be several factors that could cause an increased count, such as not cleaning your clean water bottle with BC-San, or not correctly following the Alpron protocol. If you are concerned about rising count, please email clinical support and we can discuss possible options moving forward



BC-SAN X 12 CODE -

3199-12
ORDER ONLINE
AICONLINE.CO.UK

Audit fail conditions > Failed Samplers

The majority of sampler tests will display bacteria growth as the "Dots" described on the previous page. However, there

are several results that are also fails that do not present themselves as dots. These conditions are described below.

Audit fail condition > Total overgrowth

Total overgrowth means the slide has been completely overgrown with bacteria and is a fail result. Usually this presents itself as a "slime" across the paddle.

In some cases, such as pic 1 to the right, the overgrowth and slime is obvious as it appears shiny and slimy under the light. Other times the slime is less obvious in natural light, but becomes apparent under a magnifying lamp. As such we always recommend checking samplers under a magnifying lamp.

Another way to check for slime is to run a tongue depressor over part of the paddle. As seen in pic 2 this will remove part of the slime and reveal total overgrowth.

In cases of total overgrowth this should be recorded as too many to count (TMTC).



Audit failure > What to do in case of audit failure

If your Red Sampler shows over 100 dots or TMTC we would first recommend a re-test to confirm it is failing consistently.

If your retest confirms the failed result, our protocol advises/recommends that you perform a biofilm removal treatment. This can be done over the weekend for best results (BRS + Code D1-1018).

Old chairs might require several BRS treatments as the biofilm will have had a long time to become hardened into the tubing of the dental unit.

If Sampler tests are still failing after several biofilm removals (typically 3), your chair may need to have the tubing replaced, as the biofilm may have become hardened into the tubing and can no longer be removed chemically.

Call clinical support to discuss your results: 01833 630393

BRS+ CODE -

D1-1018
ORDER ONLINE
AICONLINE.CO.UK



Auditing > Photographing your Samplers

To get a more accurate view of your Sampler results, we recommend taking photographs. Photographs under an illuminated magnifying lamp can make it easier to identify colonies and overgrowth.

If you are unsure of your results, you can also send your photos to csd@aiconline.co.uk and our team can cast a helping eye over your results. You may also be asked for Sampler photos as part of our validation process. Please follow the below method for photographing Samplers.

Step 1 Position under magnifying lamp

Take your photographs under a illuminated magnifying lamp, with the outer plastic casing off. This helps us see any colonies on the paddle. It may be helpful to have someone hold the Sampler under the lamp while you take the photo or visa versa.

Step 2 Minimise light reflection

Try to minimise the light reflection from the mag lamp. Sometimes reflection from the mag lamp can give the appearance of slime on the Sampler, when in actual fact it is clear.

Step 3 Focus

Make sure the Sampler is in focus and not blurred. We need to be able to see and count any colonies.


Step 4 Upload and rename

Rename your photo files with the name of the water source (e.g surgery 1) or provide a way we can identify the source. This is so we can determine any actions required, or if the Samplers are clear, issue your certificate!

Step 5 Assessing photos

Once renamed you can now view your photos on your computer. Use the zoom function on your desktop photo viewing software to get an excellent view of the Sampler and colonies.

Pic 1 Easy to read Sampler photo

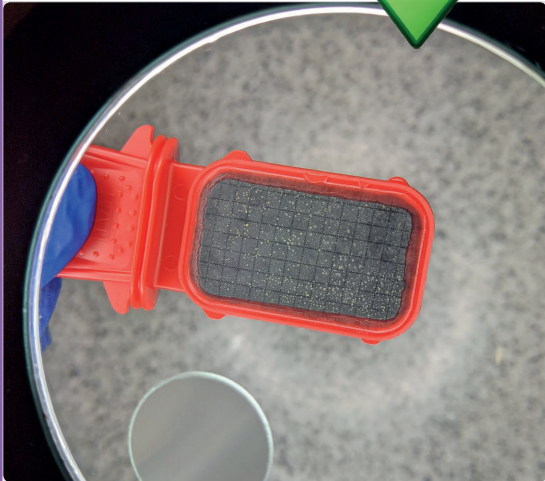

Pic 2 Difficult to read Sampler photo



Auditing > Sending your photos to Aura Infection Control





If you are unsure of your Sampler results or would like some advice, email photos and Sampler results to csd@aiconline.co.uk.

Our clinical support team will help you interpret results, advise on actions to take and issue your certificates.

Sampler Reading > Gallery of results

We see many types of Sampler results. Colonies come in many different shapes, sizes and colours, however these factors are not important by themselves, only total number of colonies.

In some cases, a slimy substance can form on the Samplers. This counts as a fail. These following images should give you some idea of what can form on Samplers to help you interpret your own results.

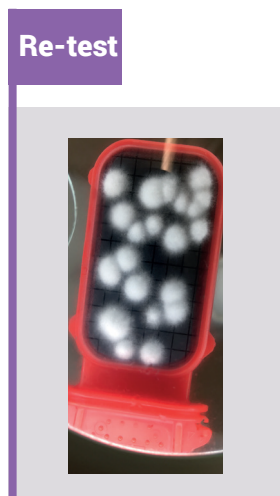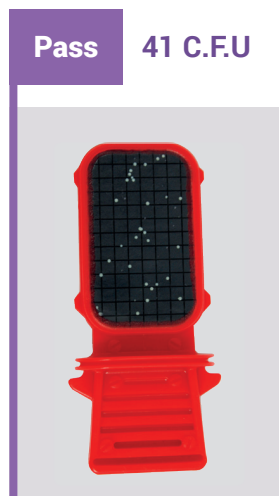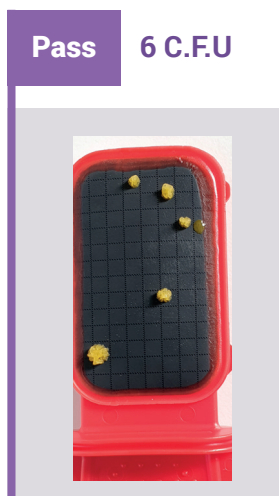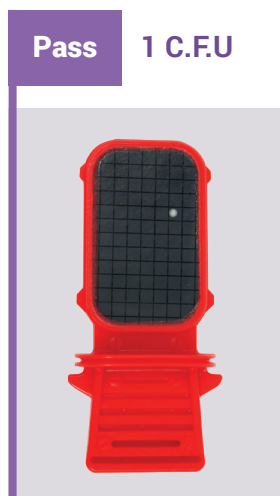

Ongoing testing > Certification and Compliance Club

Once your Sampler tests are passing, you are now ready to start routine monitoring of your dental unit waterlines. We recommend you perform a Sampler test every 3 months. Regular testing provides proof your dental units are producing drinking quality water, and alert you to any potential biofilm growth.

Our Compliance Club makes routine monitoring easy. Samplers arrive every 3 months, in time for testing, without you having to worry about ordering. Find out more about Compliance Club on page 15

We also offer certification for dental chairs that pass Sampler tests. Your water quality certificate provides evidence of proactive dental unit waterline management.

To apply for your certificate, complete our audit sheet in full (example on the next page) and return to **csd@aiconline.co.uk**.



Example audit sheet > How to complete the sheet

Once you have completed your tests and have Samplers showing less than 100 C.F.U, we can issue a water quality certificate for your records.
This is a quarterly certificate. To ensure continued certification

we will require the results of quarterly Sampler testing and **proof of purchase of Alpron and Samplers** (except where Samplers are purchased from us). Use this audit booklet to keep track of your results.

Please provide practice name, postcode and surgery name. This helps our team input the correct data for your certificate.

Sampler Quarterly Testing			
Practice Name: <i>Smith Dental Surgery</i>			
Post Code: <i>DF8 9KL</i>			
Surgery Name: <i>Surgery 1</i>			
Q1: <i>January</i>	Q2: <i>April</i>	Q3:	Q4:
Alpron Lot / Batch No. <i>501234</i>	Alpron Lot / Batch No. <i>501234</i>	Alpron Lot / Batch No.	
Sampler test date <i>13/01/22</i>	Sampler test date <i>08/04/22</i>	Sampler test date	
Sampler read date <i>20/01/22</i>	Sampler read date <i>15/04/22</i>	Sampler read date	
C.F.U. Count <i>10</i>	C.F.U. Count <i>105</i>	C.F.U. Count	
PASS	FAIL	PASS	FAIL
☒	☐	☐	☒
Required: sign your Audit Data. By signing, you are confirming information is collected following Aura's protocol and is correct.		Required: sign your Audit Data. By signing, you are confirming information is collected following Aura's protocol and is correct.	
<i>MR SMITH</i>		<i>MR SMITH</i>	

The 6 digit no. starting with a 5 is located on the side or bottom of the bottle. This acts as proof of purchase which we need for your certificate.

You should read your samplers on the 7th day of incubation i.e. if taken on 1st of the month, they should be read on the 8th.

We require the exact number of C.F.U for certification (i.e a whole number). Please count the dots on the paddle. For total overgrowth, write TMTC (too many to count).

0 - 80 C.F.U is a pass. 80 - 99 is a caution but still a pass. 100+ is a fail.

Important - Declaration
All audit data must be signed off by a responsible person. By signing, you are confirming all data is correct and collected following Aura's protocol. False information invalidates certificates.

Example audit sheet > Completing fail actions and re-tests

In the event of a Sampler fail, please complete this section, if not proceed straight to declaration. If your tests have failed or you are performing a re-test you will need to complete the fail actions section.

Corrective action taken in response to a fail, for example a BRS or BRS + Kit, or if you are re-testing only.

If you have performed a Biofilm Removal - Re-tests following a corrective BRS treatment should be taken 7 days after treatment completed. If the re-test passes, continue to declaration. If another fail is recorded, please keep this sheet and contact us for next steps

If you are re-testing only - You can perform the Sampler test as soon as you receive it.
Once again, read the Sampler after 7 days and record the exact C.F.U

Record the re-test result. If the re-test fails after a BRS please contact us for further guidance.

Once completed, please sign off your audit data.

Failed Quarter	Action taken (BRS, BRS +, Sampler re-test)
<i>April</i>	<i>BRS and Re-test</i>
Sampler re-test date	<i>11/05/22</i>
Sampler read date	<i>18/05/22</i>
C.F.U. Count	<i>20</i>
PASS	FAIL
☒	☐
Required: sign off re-test data By signing, you are confirming information is collected following Aura's protocol and is correct.	<i>MR Smith</i>

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS	FAIL	PASS	FAIL	PASS	FAIL	PASS	FAIL
☐	☐	☐	☐	☐	☐	☐	☐
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (BRS, BRS +, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (BRS, BRS +, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS	FAIL	PASS	FAIL	PASS	FAIL	PASS	FAIL
☐	☐	☐	☐	☐	☐	☐	☐
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal,, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (BRS, BRS +, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal,, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (BRS, BRS +, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Sampler Audit Sheet

Sampler Quarterly Testing							
Practice Name:							
Post Code							
Surgery Name:							
Q1:		Q2:		Q3:		Q4:	
Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.		Alpron Lot / Batch No.	
Sampler test date		Sampler test date		Sampler test date		Sampler test date	
Sampler read date		Sampler read date		Sampler read date		Sampler read date	
C.F.U. Count		C.F.U. Count		C.F.U. Count		C.F.U. Count	
PASS		FAIL		PASS		FAIL	
☐		☐		☐		☐	
Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.		Required: sign off this months Audit Data.	
By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.		By signing, you are confirming information is collected following Aura's protocol and is correct.	

Audit Sheet > Fail actions and re-tests

Only complete the following section if you have a fail result (**over 100 c.f.u.**) or you are performing a Sampler re-test. Sampler re-tests should be taken **7 days** after corrective BRS treatment. If the re-test passes, please sign off your audit data. If another fail is recorded, please keep this sheet contact us for next steps.

If no corrective actions are needed, please ignore this section.

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	

Failed Quarter	Action Taken (Biofilm Removal, Sampler re-test)
Sampler re-test date	
Sampler read date	
C.F.U. Count	
PASS	FAIL
☐	☐
Required: sign off re-test data	
By signing, you are confirming information is collected following Aura's protocol and is correct.	



DUWL Auditing Made Easy!

Simplify your audit routine with Compliance Club

The Compliance Club is our FREE service that assists practices in maintaining a quarterly audit routine for dental unit waterlines. Quarterly auditing shows you are being proactive with your dental unit waterline management, and provides evidence your dental units are producing drinking quality water.

You pay for Samplers, BC-San and Alpron at a reduced rate, saving money and reducing human error (I forgot to order them).

Benefits

- Samplers arrive automatically
- Eliminates Human Error ("I Forgot to Order Them")
- Peace of Mind
- QWS Water Quality Certificate to display in practice
- Dedicated Clinical Support Team
- 15% Discount off the full infection control range
- Certification Fee & Subscription fee: £0.00

// *Aura infection control are a great company to use for water quality control. I run a dental surgery and they have been so helpful in helping me maintain my dental unit waterlines and helping me ensure that my patient and staff safety is taken care of." T. Martin* //

Membership Levels

You can choose to add as many Samplers as you need onto your Compliance Club package. In addition you can also add other DUWL essentials on.



CC Standard

1 Sampler per chair, per quarter. You can also add additional samplers for additional test points, such as RO machines.



CC Plus

1 Sampler and 12 BC-San bottle cleaner per chair, per quarter.



CC Gold

1 Sampler 12 BC-San, and 1 bottle of Alpron per chair per quarter.

Price Comparison

Price Comparison - 1 Chair Practice

Non Member	CC Gold Member
1 x Sampler	1x Sampler
12 x BC-SAN	12 x BC-SAN
1 x Alpron	1 x Alpron
Total Cost = £138.45	Total Cost = £85.43*

To join or find out more please email csd@aiconline.co.uk



Aura Infection Control Ltd

14b Redwell Court
Harmire Enterprise Park
Barnard Castle
County Durham
DL12 8BN

Tel: 01833 630393

Web: www.aiconline.co.uk

VAT Reg: GB 181 4500 29
Company Number: 08814266